D

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Vertrieb

JW Sales GmbH

Brand Division Cosmedico Medical System

Kölner Strasse 8

D-70376 Stuttgart

Tel. +49 (0)711 54 00 4 - 0


Fax +49 (0)711 54 00 4 - 55

info@cosmedicomed.de

www.cosmedicomed.de

Instruction manual

Attention: Read these instructions carefully before using the appliance! Keep this manual for future reference.

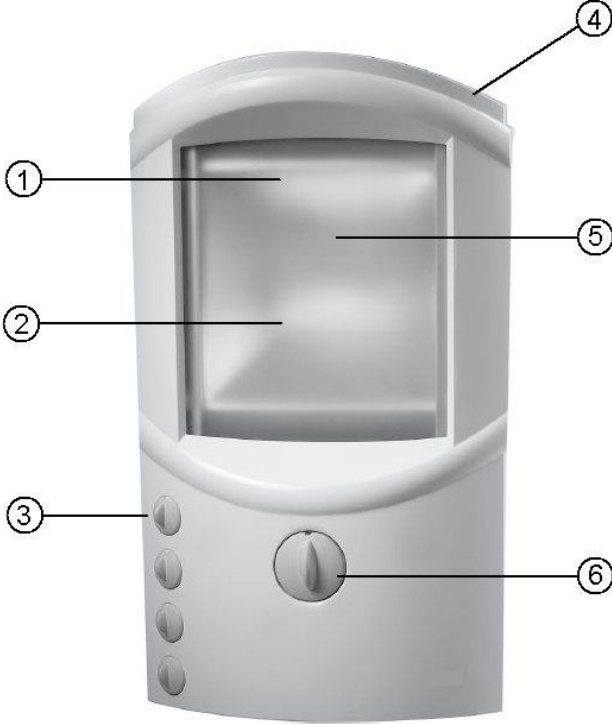
1. DELIVERY CONTENTS

COLLAGEN Pro Beauty Appliance exists of:

- 1 Filter screen
- 2 Reflector
- 3 Memo buttons for 4 users
- 4 Carrying handle
- 5 Collagen high-pressure lamp
- 6 30 minute timer

Accessories

- 7 1 pair of safety goggles
- 8 HYALURON Pro Beauty Fluid 30 ml

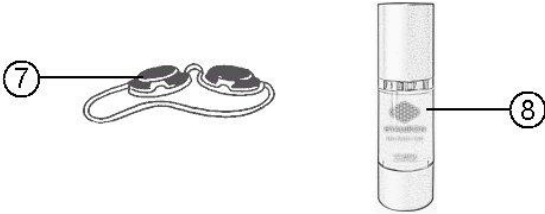2. IMPORTANT NOTES

1. Read these instructions carefully before using your device. Always follow the safety rules and instructions.
2. **Note:** None-users, especially children, should not be present when the appliance is in use.
People who are not - due to their physical, sensory or mental capabilities or their inexperience or ignorance - in a position to use the equipment safely should not use this device without supervision or instruction by a responsible person. Children should be supervised to ensure that they do not play with the device.
3. Supply only 230-volt alternating current to the appliance.
4. Do not use the appliance with a defective timer.
5. For reasons of electric safety, do not use water on the appliance or immerse it in water. Do not use the appliance in damp rooms or near a sink.
6. Caution, the lamp will not illuminate when it is hot. Let it cool down 5 to 6 minutes before using it again.
7. The lamp has a lifetime of approx. 500 hours. A lamp exchange is necessary after approx. 2700 exposures. Should you need a new lamp, please contact the customer service.
8. Make sure the filter screen does not show any signs of damage. A defective filter screen could cause serious skin burns.

9. Do not use the appliance in an explosion-hazard area.
10. Do not cover the vent slots at the back of the appliance while in use. Make sure that the vent slots are not blocked by impurities.
11. Unplug the appliance after use.
12. Only store the appliance after it has cooled down (cooling period ca. 15 minutes).
13. If a defect of the appliance should occur, please contact the customer service.

3. WHAT YOU SHOULD KNOW ABOUT COLLAGEN AND HYALURONAN

With light wave appliance COLLAGEN Pro Beauty and HYALURON Pro Beauty Fluid by Cosmedico Medical Systems you can rejuvenate your skin visibly and perceptibly. Studies have shown that the body's own collagen production can be stimulated with red light of a special wavelength range. The COLLAGEN Pro Beauty Appliance emits red light waves of the ideal wavelength range. The so-called Photo Bio Stimulation in the deeper layers of skin activates the metabolism of the skin and leads to a reactivation of the collagen synthesis in the skin. The ageing skin is consequently stimulated to produce skin-lifting collagen and elastin again.

The combination of the red light waves with the HYALURON Pro Beauty Fluid ensures an optimal action endorsement. Before each light wave treatment, the fluid is amply put on the face. The light wave exposure channels the active ingredients contained in the fluid into the deeper layers of skin. This elevates the level of hyaluronic acid in the skin to a teenage level. The elasticity of the connective tissue is heightened, the skin stability rejuvenates. At the same time the collagen synthesis of the skin is stimulated.

The use of the HYALURON Pro Beauty Fluids in combination with the COLLAGEN Pro Beauty Light Wave Appliance clearly strengthens the rejuvenating effect and leads to a visible smoothing effect on the skin surface. The result is a firmer and tighter skin with reduced wrinkles and a fresh, young look.

APPLICATION ADVANTAGES:

- Makes skin rejuvenation possible at home with visible anti-aging effect
- Reduction of wrinkles (face, décollete)
- Improvement of the skin elasticity, suppleness and resilience
- Stimulation of the collagen synthesis of the skin
- Activation of the skin metabolism
- Boost of a fresh, young look of the skin
- Heightening of wellness, relaxation effect for body and soul
- No side effects

4. INSTALLATION OF THE APPLIANCE

COLLAGEN Pro Beauty has been developed as a table-top device. You can put the appliance on a suitable table. Your face should be at the same height as the red filter screen. The distance between face and filter screen should be at least **15 cm**.

5. USE

The use of COLLAGEN Pro Beauty is recommended especially for people with mature skin. The following intervals have proved to be very effective for therapy with red light waves.

FREQUENCY OF USE AND DURATION

Exposure distance: **15 cm** to the red filter screen of the appliance.

Keep your eyes closed during use.

 **Always use the goggles to protect your eyes during exposure.**

If you need another pair of safety goggles, please contact the manufacturer and customer service.

Remove the make-up of your face before each treatment and clean your face, and put the HYALURON Pro Beauty Fluid evenly on the skin.

GUIDE OF USE		
Exposure distance 15 cm		
Week	Number of treatments	Duration per treatment
Week 1	2	12 min.
Week 2	2	12 min.
Week 3	2	12 min.
Week 4	2	12 min.
Week 5	2	12 min.
Week 6	2	12 min.
Week 7	2	12 min.
Week 8	2	12 min.
Week 9	3	12 min.
Week 10	3	12 min.
Week 11	3	12 min.
Week 12	3	12 min.
Week 13	0 (Pause)	0
Week 14	0 (Pause)	0
Week 15	0 (Pause)	0
Week 16	0 (Pause)	0
Continuation see below		
Max. number of treatments per year: 90		

When the month-long pause has ended, you can continue the treatment 2 times per week for 2 months. After the 2 months, pause again for 1 month. And so on. The pauses have as purpose not to weaken the skin in its own function, but to prompt it to self-activation and regeneration. In case of a continuous exposure, the own mechanism of action of the skin could decrease.

A yearly number of 90 exposures in total should not be exceeded.

6. STARTING THE APPLIANCE

Plug in the table-top device. Turn the timer knob clockwise to set the exposure time. This starts the appliance and the timer. The system stops automatically at the end of the exposure time.

7. EXPOSURE ADVICE

Keep an exposure distance of 15 cm. The face should be centered at the height of the red filter screen. Be sure to put the goggles on. If you like, you can keep your eyes closed under the goggles. A light alternating motion of the head does not cause any disadvantage. The decollete can be treated as well.

8. MEMO BUTTONS

The number of already performed exposures can be tracked and checked if you set the memo buttons. Four memo buttons for up to four users are available. For each use, the corresponding user turns his memo button one step to the right.

9. CLEANING AND CARE

Pull the mains lead before starting to clean. Always make sure the filter screen remains clean and intact. Spots, dust and fingerprints can burn and reduce the efficiency of the radiations. To clean the screen, use a clean cloth and glass cleaner (never use white spirit or household gasoline). If necessary, clean the housing with a soft damp (but not wet) sponge or cloth and a mild detergent.

10. REPAIRS

- The filter screen has special filtration features and in case of damage, it should be replaced by an original filter screen.
- Broken glass is excluded from warranty, even on new appliances. Likewise, no warranty can be claimed for damage in transit or repaired systems that are poorly packed.
- We advise you to keep the original packaging materials.
- All repairs, including the replacement of the power cord, should be made by the manufacturer or his customer service or a qualified technician, in order to avoid danger.

PROTECTION OF THE ENVIRONMENT – DIRECTIVE 2002/96/EC

In order to preserve our environment and protect human health, the waste electrical and electronic equipment should be disposed of in accordance with specific rules with the implication of both suppliers and users. For this reason, as indi-

cated by the symbol  on the rating label or on the packaging, your appliance should not be disposed of as unsorted municipal waste. The user has the right to bring it free of charge to a municipal collection point performing waste recovery by means of reuse, recycling or use for other applications in accordance with the directive.

TERMS OF GUARANTEE

For our products, the legal terms concerning claims under guarantee are valid. During the legal period, all material and manufacturing defects are covered. No claims will be accepted for broken accessories, defects resulting from inappropriate use or care, excessive wear due to intensive use, wrong use due to the non-observance of the instruction manual,

or normal wear and tear. No guarantee claims will be accepted if the appliance has been modified or if attempts have been made to repair the appliance. In such case, the vendor cannot be held responsible in case of accidents.

In the event of claim, send the appliance cleaned, well packed and with a short, easily readable description of the default as well as the receipt freight prepaid to the customer service. Claims under guarantee will only be valid if the receipt or a copy of a repair bill is handed over with the appliance that needs repair.

AFTER THE LEGAL PERIOD OF GUARANTEE

Repairs after the guarantee period can be performed at a charge by a competent qualified technician or a repair service.

Manufacturer and customer service

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Gustav-Töpfer-Straße 6

07422 Bad Blankenburg / Thüringen

Germany

Phone Germany 036741 70268 International 0049 36741 70268

Fax Germany 036741 70227 International 0049 36741 70227

service@efbe-schott.de

www.efbe-schott.de

Distributor

JW Sales GmbH

Brand Division Cosmedico Medical System

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Efbe Elektrogeräte GmbH